

## Recombinant Human Interleukin 18 Binding Protein (IL18BP)

Catalog No.: TD-TP02620 50µg

### Sequence Information

**Species:** Human

**Gene ID:**10068

**Swiss Prot:**O95998

**Synonyms:**Interleukin-18-binding protein,  
IL18BPα

**Residues:**Thr31~Gly194

TPVSQTTTAAATASVRSTKDCPCSPQPPVFPAAKQCPALEVTWPEVEVPLNGTSL  
SCVACSRFPNFSILYWLGNFSFIEHLPGRWEGSTSRERGSTGTQLCKALVLEQ  
LTPALHSTNFCVLDPEQVVQRHVLAQLWAGLRATLPPTQEALPSSHSSPQQ  
QG

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Secreted

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.5

**Predicted Molecular Mass:** 18.9kDa

**Accurate Molecular Mass:** 19kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

电话: 027-87561633 手机: 189-8623-2100

QQ客服: 121545059 网址:<http://www.itiande-neuron.com>

公司地址: 武汉东湖新技术开发区高新大道666号光谷生物城生物创新园B1栋b274-9室

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

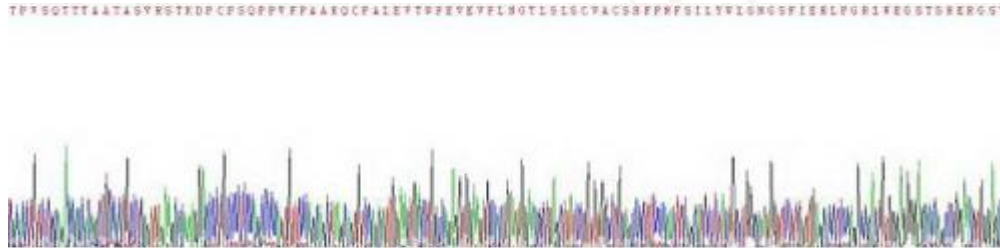


Figure 1. Gene Sequencing (Extract)

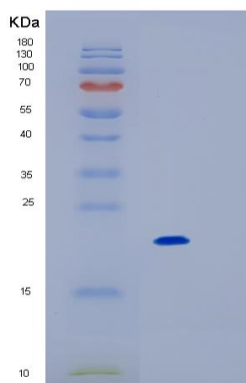


Figure 2. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.