

## Recombinant Cluster Of Differentiation 58 (CD58)

Catalog No.: TD-TP01523 50µg

### Sequence Information

**Species:** Human

**Gene ID:**965

**Swiss Prot:**P19256

**Synonyms:**LFA3; Ag3; Surface glycoprotein

LFA-3; Lymphocyte Function Associated

Antigen 3

**Residues:**Met1~His214

MVAGSDAGRALGVLSVVCLLHCFGFISCFSSQIYGVVYGNVTFHVPSNVPLKEV

LWKKQKDKVAELENSEFRAFSSFKNRVYLDTVSGSLTIYNLTSSDEDEYEMESP

NITDTMKFFLYVLESLSPTLTCALTNGSIEVQCMPIPEHYNSHRGLIMYSWDCP

MEQCKRNSTSIYFKMENDLPQKIQCTLSNPLFNTTSSIIILTTCCIPSSGHSRH

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** Two N-terminal Tags, His-tag and GST-tag

**Subcellular Location:** Nucleus.

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.7

**Predicted Molecular Mass:** 54.1kDa

**Accurate Molecular Mass:** 54kDa as determined by SDS-PAGE reducing conditions.

### **[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### **[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

电话: 027-87561633 手机: 189-8623-2100

QQ客服: 121545059 网址: <http://www.itiande-neuron.com>

公司地址: 武汉东湖新技术开发区高新大道666号光谷生物城生物创新园B1栋b274-9室

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### [ IDENTIFICATION ]

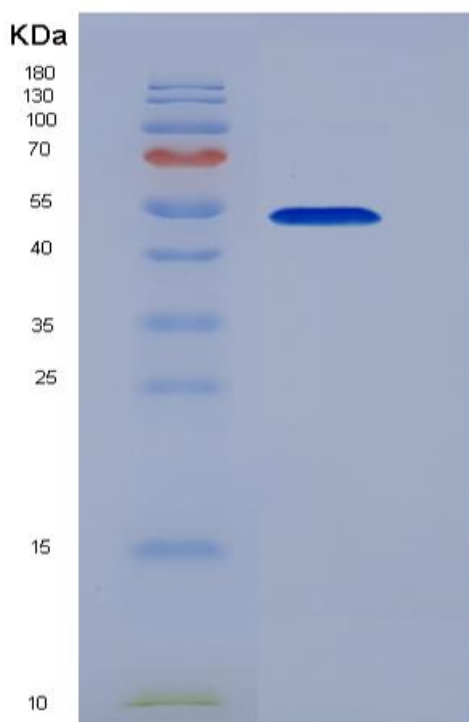


Figure 1. SDS-PAGE

### [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.